

Biomedical



Improving  
efficiency  
in cell culture  
and upscaling

Dissolvable  
microcarriers


[rousselot.com](https://rousselot.com)

**Rousselot** |  **IamFluidics**  
| by Darling Ingredients revolutionize microparticles

As demand for cell culturing systems continues to rise — driven by applications in personalized medicine, cell therapies, biologics, and cultivated meat — there is a growing need for efficient and scalable cell culture systems.

## The new dissolvable microcarrier makes scalable cell culturing efficient


The dissolvable microcarrier (DMC), developed through the collaboration between Rousselot® and lamFluidics®, unlocks the full potential of cell culturing in bioreactors. By providing an unmatched efficiency in cell recovery paired with seamless integration into standard cell culture processes, our dissolvable microcarrier will drastically improve process efficiency.

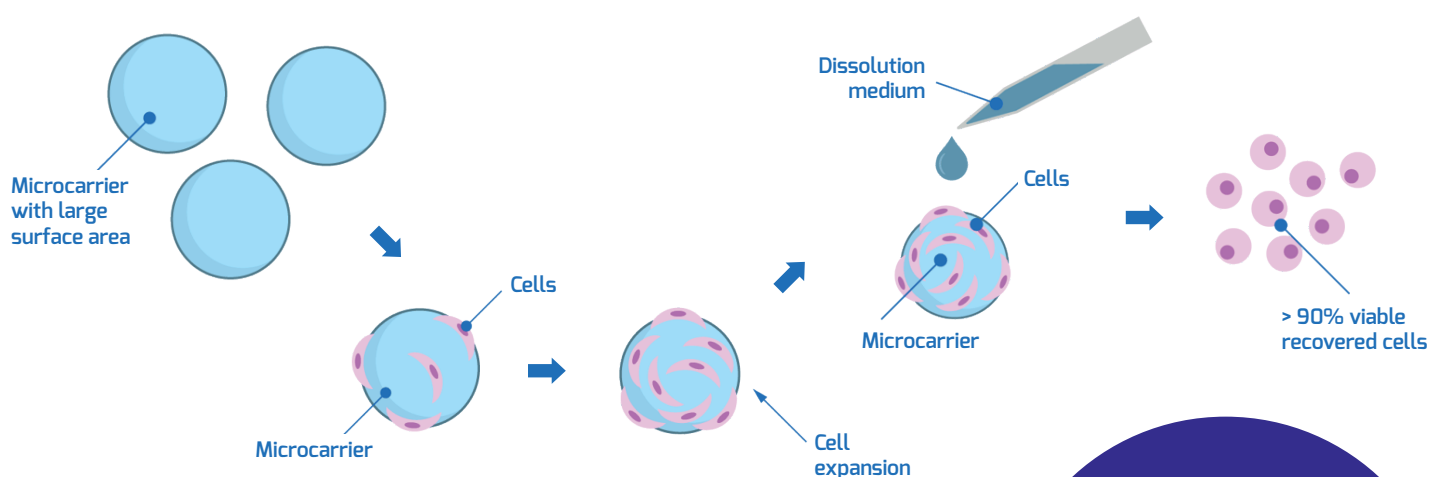

### Features

- + **Dissolvable** within 15 minutes
- + **High harvesting yield** (> 90%)
- + **Transparent for easy visualization** of cells during culture
- + **Collagen (denatured) coating** – purity controlled and excellent batch-to-batch consistency
- + **Full traceability** and documentation (virus inactivation) available
- + **Easy integration into standard cell culture procedures**
- + **Uniform size & shape** (CV < 10%) for consistent sedimentation & seeding
- + **Physiological soft cell-material interface**
- + **100% free of microplastics**
- + **Compatible with widely used adherent cell types** (e.g., Mesenchymal stromal cells, iPSC/ESC-derived cells, and producer cell lines like VERO, CHO & HEK)


Harnessing the biocompatibility of naturally derived collagen while safeguarding against animal-derived pathogen transmission

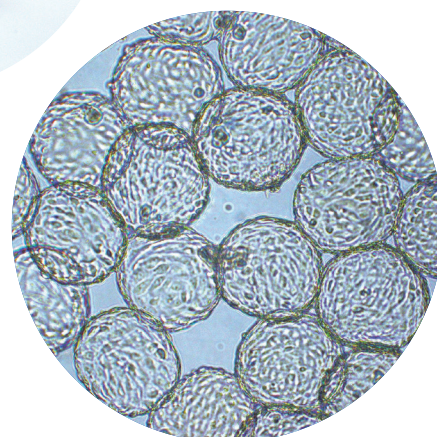



## Enhanced cell attachment and proliferation



The denatured collagen coating promotes cell attachment and proliferation, and supports physiological interactions between cells and the material. The dissolvable microcarriers are suitable for mesenchymal stem cells (MSCs) and other adherent cell lines.

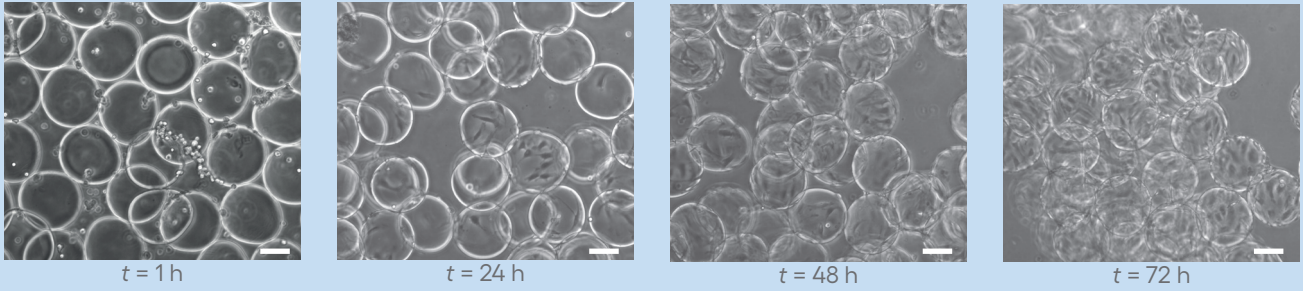

hMSCs (human mesenchymal stem cells) grown on dissolvable microcarriers show uniform cell attachment within the first 24 hours and confluency after 72 hours. This proves that the MSC's have a great bead-to-bead transfer on these DMC's.



## Readily dissolvable using standard culture reagents

The microcarriers can dissolve within 15 minutes after the addition of standard cell culturing reagents such as EDTA/TrypLE™ or Trypsin, allowing for efficient cell harvesting, cell recovery and swift integration into existing cell-expansion systems.

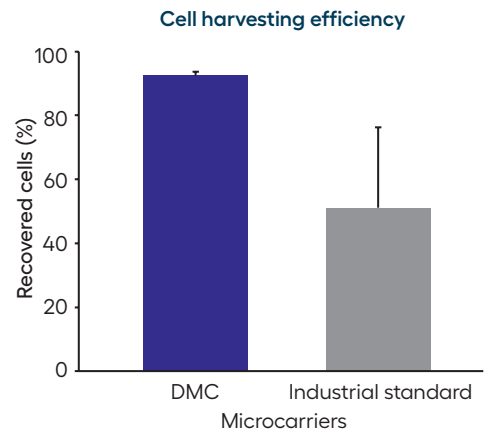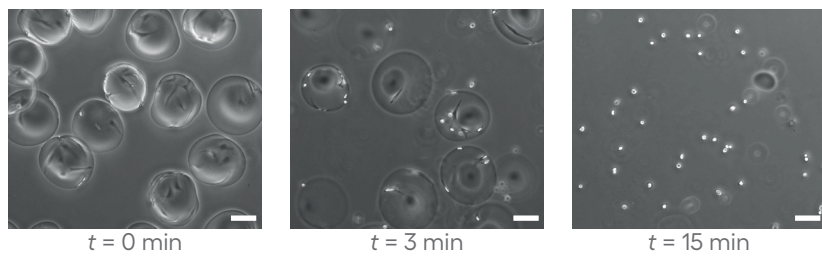

## Over 90% cell viability

With over 90% cell harvesting efficiency and viability, the dissolvable microcarriers dramatically reduce cell loss compared to the typical 50% seen with non-dissolvable microcarriers.

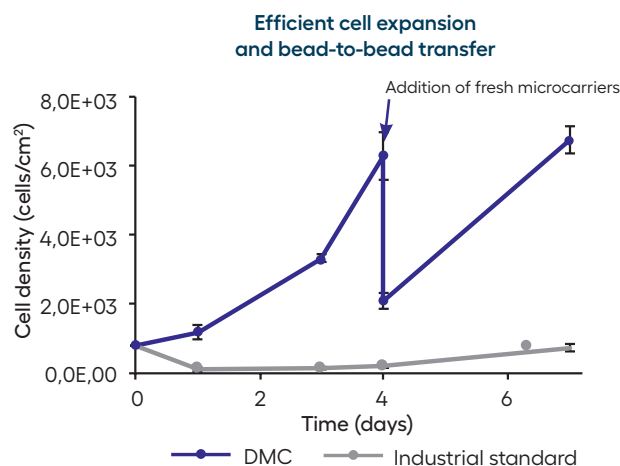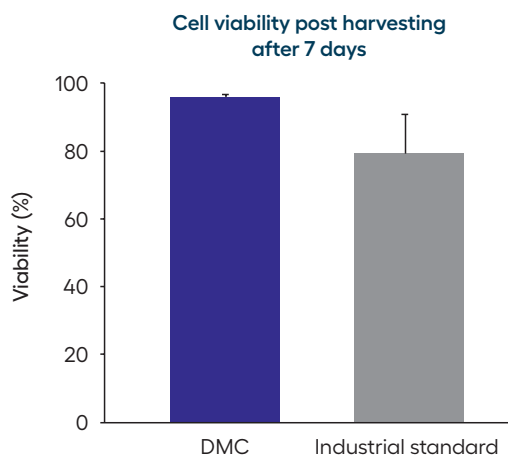

In a bioreactor, cells expanded for 7 days with surface area expansion at day 4, showed excellent cell viability above 97%.

## Monodisperse & transparent for consistent and easy handling

The monodisperse and transparent microcarriers enable reproducible and easy handling, following standard cell culturing protocols.

## Product range

We offer the dissolvable microcarriers as dry powder standard in 1 gram and 10 grams.

Larger quantities (100 grams and 1 kg) are available on request.

Available online  
[www.microparticles.shop](http://www.microparticles.shop)



Scan QR code for more information, cooperation and questions on GMP.

Rousselot can help you with virtually any product requirement or innovation you have in mind, offering:



Transparency



Full traceability



High standards of quality and safety



Committed to the environment and to our clients





Global support and expert advice

### Your sales contact information



Rousselot, Darling Ingredients' premier collagen and gelatin brand, and lamFluidics, a microparticle innovator, have partnered to launch dissolvable microcarriers for adherent cell culturing. Combining Rousselot's premium collagen expertise with lamFluidics' patented IN-AIR MICROFLUIDICS™ technology, the microcarriers enhance cell attachment, rapidly dissolve for efficient harvesting, and offer over 90% cell viability. This scalable, cost-effective solution addresses industry challenges in advanced cell therapies, biologics manufacturing, and cultured meat production, positioning both companies as leaders in biopharmaceutical innovation.

 Rousselot  
 @RousselotBio

**Rousselot B.V.**  
Kanaaldijk Noord 20  
5681 NM Son  
The Netherlands

+31 499 364 100  
[rousselot.com](http://rousselot.com)

**Rousselot**  
| by Darling Ingredients

**Disclaimer** All rights reserved. No part of this brochure may be reproduced, distributed or translated in any form or by any means, or stored in a database or retrieval system, without the prior written permission of Rousselot. Rousselot alone retains the copyright to the entire content of this brochure and the intellectual property rights to all designations of our products stated in this brochure and intellectual property rights to the products themselves. Nothing in this brochure constitutes a license (explicit or implicit) of any of Rousselot's intellectual property rights. The duplication or use of product designations, images, graphics and texts is not permitted without Rousselot's explicit prior written consent. Rousselot makes no representation or warranty, whether expressed or implied, of the accuracy, reliability, or completeness of the information, nor does it assume any legal liability, whether direct or indirect, of any information. Use of this information shall be at your discretion and risk. Nothing herein relieves you from carrying out your own suitability determinations and tests and from your obligation to comply with all applicable laws and regulations and to observe all third party rights. This product is not intended to diagnose, treat, cure, or prevent any disease. You should always consult your medical provider when using the product together with medical treatments, diets or fitness programs. The uses and claims for Rousselot's products recommended in the brochure should be adapted to the current local regulatory environment. This statement has not been evaluated by the Food and Drug Administration.